*Photos are non-contractual - they may differ from actual products

NAB TECH

FLASH POINT TESTER ABEL CLOSED-CUP METHOD



STANDARDS

ISO 13736, IP 170 and related methods.

SCOPE

These test methods cover procedures for determination of the manual and automated closed-cup flash point of combustible liquids having flash points between $-30.0~^{\circ}\text{C}$ to $75.0~^{\circ}\text{C}$. However, the precision given for this method is only valid for flash points in the range $-8.5~^{\circ}\text{C}$ to $75.0~^{\circ}\text{C}$.





NAB Tech is a flash point tester developed by Normalab. Normalab has been designing and manufacturing petroleum testing instruments and glassware since 1963.



SPECIFICATIONS

- Temperature range from -30°C to 110°C (cryostat is mandatory for temperatures below ambient)
- Temperature display in °C or °F
- 4.3" resistive color touchscreen
- Class B Pt 100 probe
- Barometric correction with temperature compensation
- Predefined temperature range 1°C/min
- Quick cooling system
- Ignitor presentation by push button
- Automatic detection by thermocouple (no calibration required)
- Electric ignitor, adjustable from 0 to 30 W (no gas)
- Sample cooling integrated
- Adjustable and settable strirring speed from 0 to 30 Rpm
- Specific treated cover against freezing (easy to clean)



CONNECTION


- 1 RS232C port
- 2 Micro USB port for update

*Photos are non-contractual - they may differ from actual products.

NABTECH V0.9 - V0.9 20.5 °C READY TO GO Bath temperature: 21.9 °C Sample name: Test DEMO Operator name: Normalab Method name: ABEL Time test: 00:00:00 Stirring speed: 0.0 rpm Expected temperature: 50.0 °C

Method name : ABEL Standard name : IP 170 Created/Modified by : Normalab Creation/Modification date : 19/07/2019

Archives Archive: 0/5 Date and time of test: - 47:47:47 Method name: ABEL Standard name: IP 170 Sample name: Operator name: Time test: 00:00:00 Barometric Pressure: 0.00 mbar Test status: Manual stop test Flash point: XXX.X °C Corrected flash point: XXXX.X °C

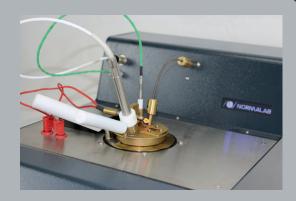

APPLICATION

- Petroleum industry
- Cosmetics and perfumes
- Pharmaceutical industry
- Chemical industry
- REACH applications
- Transportation class for hazardous goods

Biodiesel

Lubricants

SOFTWARE FEATURES

- User-friendly and intuitive software
- Predefined file
- Predefined temperature ramp (1°C/min)
- Simplified temperature calibration
- Up to 5 calibration offsets
- Software update through micro USB
- Multilingual (FR / EN)
- Data transfer or printer by RS 232
- 20 last results in memory

SAFETY FEATURES

- Access protected by 1 password
- Protection of the electric ignitor by thermal fuse
- Removable tray with drain if sample overflow or condensate at low temperature
- Dedicated alarm for each error message
- Sound alarms
- Visual indicator

OPTIONS

- Cryostat (needed below the ambient)
- RS 232 ticket printer
- NCAL, calibration kit
- Lid and cover support





Calibration

ORDERING INFORMATION

NAB Tech is a compact instrument providing quick and reliable results. Operation is made easy and safe with minimum operators intervention.

Scope of delivery:

42300 The NAB Tech tester includes

- Lid with ignitor and thermocouple
- Cup
- Stirring cable
- Power cable
- Stylus

Site requirements:

- Power supply: 230 V (+/- 10%), 50/60 Hz 4 A
 Dimension: (W) 280 x (D) 440 x (H) 250 mm
- Weight: 15 kg

SUMMARY

Temperature range	- 30 à + 110°C	Detection	Thermocouple
Temperature accuracy	0,2 °C	Ignition	100% electrical
Temperature unit	°C ou °F	Cooling system	Mecanical fan
Temperature measurement	Probe Pt 100 (Class B)	Data storage	Micro SD (20 results)



CONTACT: sales@normalab.com

Normalab FRANCE SAS ZA Caux Multipôles 1 - 76190 Valliquerville

Tel.: +33 232.700.100 Fax: +33 232.704.732 **DISTRIBUTED BY**