



Water determination in gases

aquamax KF
PRO LPG




Conform
to Standard
ASTM
D 7995 - 19

aquamax KF

PRO LPG

Description

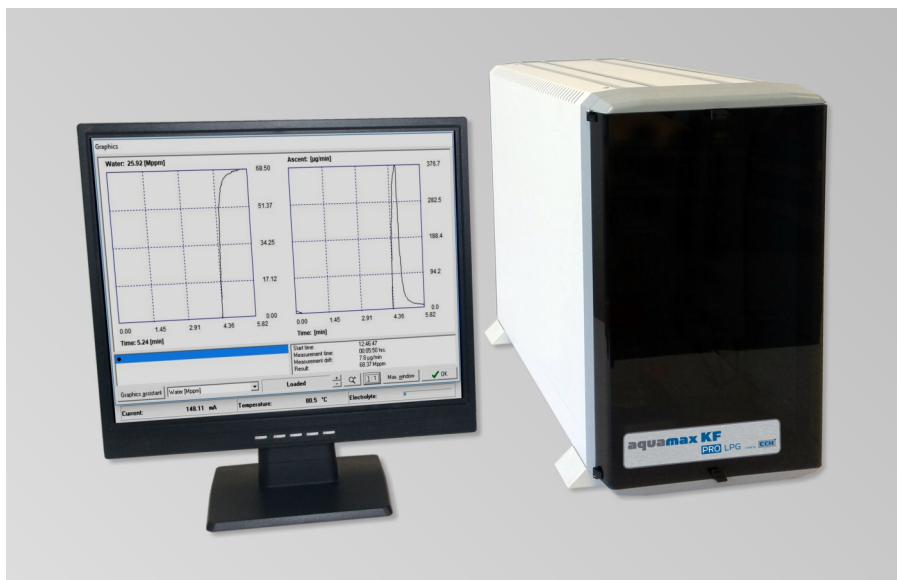

The Aquamax KF PRO LPG fulfils the requirements of the standard ASTM D 7995 - 19: Standard Test Method for Total Water in Liquid Butane by Liquefied Gas Sampler and Coulometric Karl Fischer Titration.

The Aquamax KF PRO LPG is designed for an easy and accurate determination of water in liquefied and gaseous samples such as LPG and LNG. The new device combines coulometric Karl Fischer method with an unique gas evaporation and dosing procedure.

The Aquamax KF PRO LPG includes all features required for ppm level water in LPG and Gas, including the sulphur removal cartridge eliminating all side reactions and our sample loop principle allowing you to fully automate the measurements, up to 125 per day!

All Aquamax KF PRO LPG parts are totally enclosed making this system completely safe and robust for use in the demanding petroleum industry.

The unique ECH sample loop allows you to use the instrument in your laboratory with full automation, as a portable/field use analyzer or can be integrated in to your process as an on-line system.

Applications

LPG, LNG:

- propane, propene, butane, butene, butadiene
- ethylene oxide
- chlorinated hydrocarbons, e. g. methylene chloride, ethylene chloride, vinyl chloride

Analysis of refrigerants:

- halogenated hydrocarbons

Analysis of permanent gases:

- natural gas
- technical gases
- mixtures of test gases



Features

- Determination of moisture in liquefied and gaseous samples
- Inlet pressure up to 200 bar/2900 psi
- Determination of pressure in the sample loop
- Automatic pressure regulation
- Transfer line with direct injection
- Automatic rinsing bypass and steps for rinsing
- Measuring cell without diaphragm (only one electrolyte required)
- Setting of application-specific methods
- Avoiding of side reactions by sulphur trap
- Type of result: μg , ppm (gas volume), Vppm, Mppm, Mol ppm by using the formula generator



Sulphur trap for elimination of H_2S and mercaptans

Results


Example of a measurement series with sulphur trap

Result overview:

Measurement	Sample amount	Result
1	539.282 mL	48.30 Mppm
2	539.067 mL	47.98 Mppm
3	539.282 mL	47.95 Mppm
4	538.563 mL	47.54 Mppm
5	538.555 mL	47.33 Mppm
6	538.141 mL	45.79 Mppm
7	536.514 mL	46.72 Mppm

Statistics:

Arithmetical mean:	47.37 Mppm
Standard deviation:	0.87 Mppm
Rel. standard deviation:	1.83 %

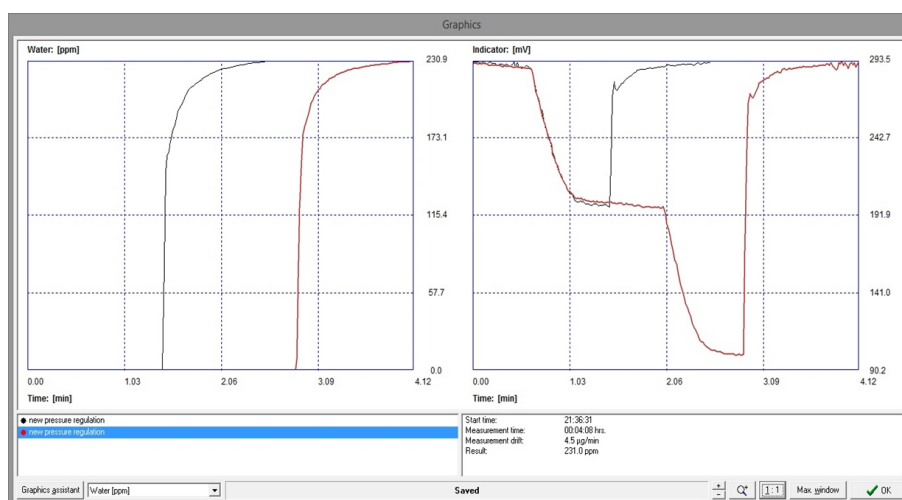

Example for multi-injection of the sample: one-step and two-step dosing process in comparison

Advantages

- Sulfur removal cartridge eliminating side reactions
- No interference calculation required
- 250 measurements can be performed in 48 hours
- Totally automated process, no operator input required for the test
- Suitable to test all gas types without any calibration or adjustments
- No separate rinsing gas is required
- Rinsing process is fully automated
- No balance is required
- High sample throughput and long reagent life
- Compact device



Typical sample cylinder with valve, e. g. of DME



Example of a 5 L Propane Cylinder

Specifications

Measurement method:	Coulometric Karl Fischer titration
Sample:	Pressurized gas sample (LNG, LPG)
Sample dosing:	Pressurized bottle or directly from the gas line
Pressure reducer:	internal (with heating element)
Sample loop:	300 mL (gas)
Rinsing and dosing:	0 ... 15 steps for each, adjustable
Measuring range:	1 µg ... 100 %
Resolution:	0.1 ppm
Detection limit:	1 ppm
Power supply:	230 V/50 Hz; 115 V/60 Hz
Dimensions:	33 x 49 x 48 cm (W x D x H)
Weight:	24 kg
Device control:	PC software (PC not included in the scope of delivery)

We are here for you



Headquarters in Germany
ECH Elektrochemie Halle GmbH
Otto-Eissfeldt-Str. 8
D-06120 Halle (Saale)
Germany
Tel.: +49 345 279570-0
Fax: +49 345 279570-99
E-mail: info@ech.de
Website: www.ech.de

Sales and Service Center in UK
ECH Scientific Limited
Building 69, Wrest Park, Silsoe
Bedfordshire, MK45 4HS
United Kingdom
Tel.: +44 1525 404747
Fax: +44 1525 404848
E-mail: info@ech.de
Website: www.ech.de