

NPM 450

AUTOMATED FLASH POINT

BY PENSKY-MARTENS CLOSED CUP TESTER



STANDARDS

ASTM D 93, IP 34, ISO 2719, NF M07019, DIN 51758, GBT 261, JIS K 2265-3 and related methods.

SCOPE

These test methods cover the determination of the flash point of petroleum products in the temperature range from 40 °C to 370 °C, and the determination of the flash point of biodiesel in the temperature range of 60 °C to 190 °C. Flash point is used in shipping and safety regulations to define flammable and combustible materials.


NPM 450 is the 5th generation of automated flash point developed by Normalab. Normalab has been designing and manufacturing petroleum testing instruments and glassware since 1963.



« EFFICIENT
AND COMPACT
UNIT »

SPECIFICATIONS

- Small size
- Capacitive touch screen (6,4 inch HD)
- One hand easy operation
- Sample cup, probe, thermocouple holder stand
- Test temperature range: 0 to 410°C
- Unique combine ignition system (electric and gas)
- Standard methods: A, B, C
- Integrated barometric correction
- Flash point detection by thermocouple
- Automatic gas flame relighting
- User programmable methods
- Heavy sample preheating and automatic stirring start
- Fast heating and cooling

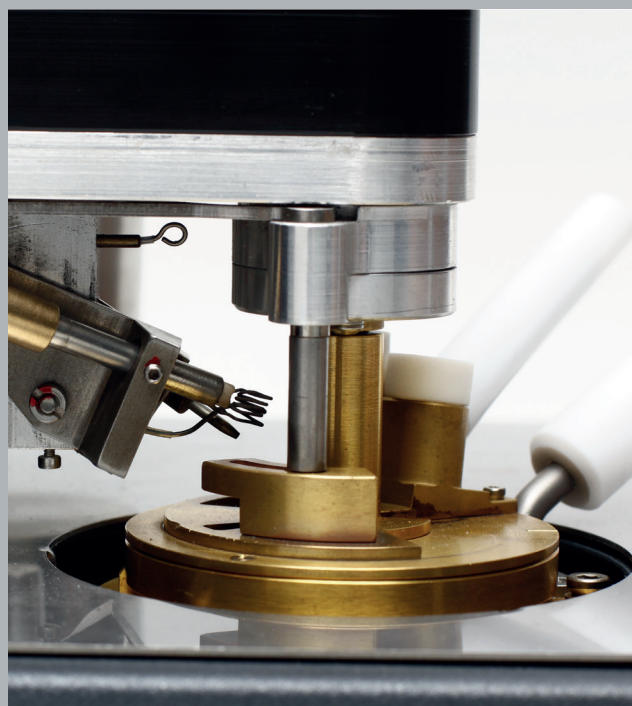

*Photos are non-contractual - they may differ from actual products.

CONNECTION

- 4 USB ports
- 1 RS232C port or LIMS
- 2 LAN ports
- 1 VGA port




SAFETY FEATURES

- Illuminated extinguisher push button
- Gas flame detection
- Over-temperature protection
- Built in gas leak detector

SOFTWARE FEATURES

- One hand easy operation
- Windows 8 operating system
- Lifetime data storage - 16 GByte SSD
- Reporting: PDF, Email, LIMS, Excel, Word...
- Different access level settable (Operator - supervisor or Maintenance...)

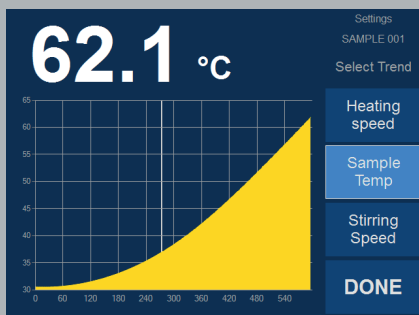

Flash Point:	0.0 °C	Test details SAMPLE 001 Select function
Corrected FP.:	0.0 °C	
Sample ID:	SAMPLE 001	Send to Trends DONE
Method:	Method A	
Exp. Flash Point:	65.0 °C	
Barometric Pressure:	992 mBar	
Min. Heating Speed:	2.9 °C/min	
Max. Heating Speed:	5.6 °C/min	
Start time:	06/04/16 08:56:57 de.	
Duration:	11 min, 14 sec	
Operator Name:	ADMIN	

EASY CONFIGURATION

Access from basic to advanced parameters made easy through the intuitive software and large screen.

- Data export through USB, LIMS, PDF or e-mails.
- User and password settings
- Data transfer to local printer
- Product files configuration

ACCESSORIES

For safety purpose, an Automatic N2 fire extinguisher option is available. This extinguishing system allows a fire to be switched off by smothering. The inert gas extinguishing system uses nitrogen. It acts on fire by reducing the oxygen content without risk for human health under normal conditions of use.

Its diffusion during emission is very fast and its density comparable to that of air (composed of 78% nitrogen) allows a durable maintenance of the extinction concentration, unlike extinguishing agents with a higher density than air.



Safety option
REF : 60410

ORDERING INFORMATION

Fully automated Pensky Martens Flash Point tester

60400

NPM 450 operable through integrated touchscreen and Windows 8 based software.

Scope of delivery:

The NPM 450 delivered ready to use with:

- Cup (P/N 60346)
- Cover (P/N 60447)
- Pt 100 probe (P/N 60453)
- Detection thermocouple cable (P/N 60448)
- Electric ignitor (P/N 60480)

Site requirements:

- Power supply: Ac 110/2400 V, 50/60 Hz - 4 A
- Dimension: (W) 260 x (D) 460 x (H) 185 mm
- Weight: 12 kg

Propane/Butane pressure: mini 20 to 40 mBar max

SUMMARY

Flash point range	0 to 410°C	Gas connections	Propane / Butane, Nitrogen
Heating rate	0 to 15°C/min	Automatic switch off	Integrated
Safety stop	0 to 30°C (above expected flash point)	Hydrocarbon detection	Integrated
Ignition	Electric and gas	Overheating protection	Integrated 600°C
Electric ignition	Regulated / monitored	Real time clock	Integrated
Stirring speeds	50 to 300 RPM (+/-3 RPM)	Fire detection	Ionization
Flash detection	Thermocouple	Extinguisher push button	Integrated & illuminated
Cooling	Integrated - Forced air	Automatic fire extinguisher	Integrated (nitrogen network or optional bottle)



CONTACT : sales@normalab.com

Normalab FRANCE SAS
ZA Caux Multipôles 1 - 76190 Valliquerville
Tel. : +33 232.700.100
Fax : +33 232.704.732

DISTRIBUTED BY